

# SPAN Advantage



- Orthopedic grid deck for flexibility and ventilation
- Low-bed design with advanced clinical positioning
- Optional tool-less expansion up to 42"

science, safety, simplicity


 savaria®



# SPAN Advantage



Easy-on expanders provide a wider surface for user comfort, fall prevention and ADLs

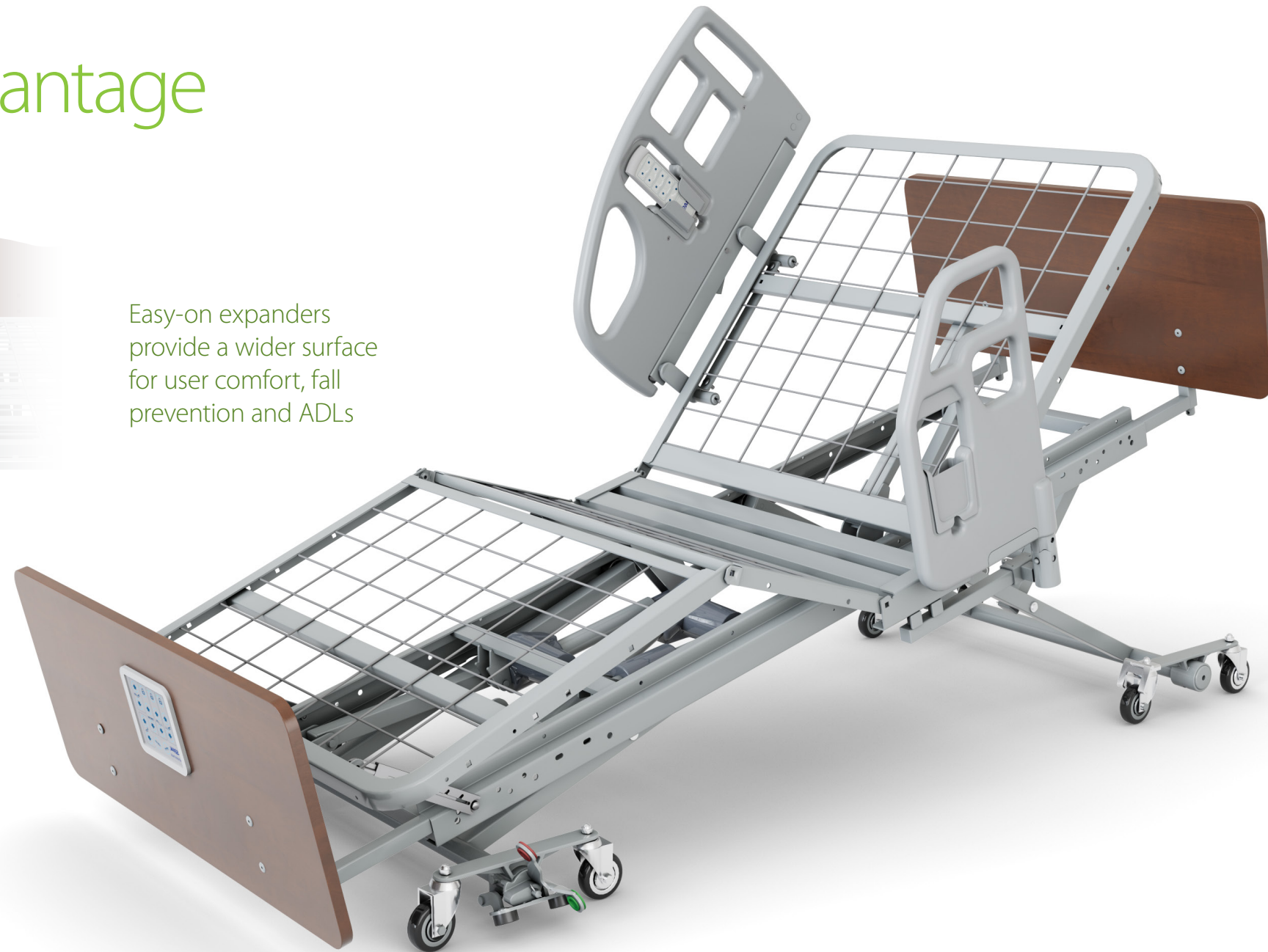

## Elevate patient care and comfort, boost caregiver safety and ergonomics

### Packed with features:


- Precision welded orthopedic grid deck for no-sag support
- Durable, antimicrobial polyester powder-coating
- Vascular footrest with easy ratchet operation
- Vertical movement with virtually no horizontal movement
- Arched legs allow wide stance for clearance and stability
- Baseboard bumpers protect walls
- Standard hand pendant controller
- Smooth deck surface, free from any protruding welds, for easy disinfection
- Anchored corner retainers prevent mattress shifting and dangerous gaps

### Customize with:

- Ergonomic 3-position side assist rail bed for sturdy ingress/egress support
- Easy, quick width adjustable, frame mounted Molded Pivot Assist pair
- Side carrying transport dollies to move beds through narrow halls and ergonomic maintenance
- Tool-less bed expanders for 42" width
- 4" bed extenders to 84" length
- Molded plastic quick release half-rail
- Foot-end controller with staff lockout

- Valley Oak  
(Standard boards only)
- Candle Light
- Hardrock Maple
- Summer Flame
- Choco

### Choice of Finish & Profile

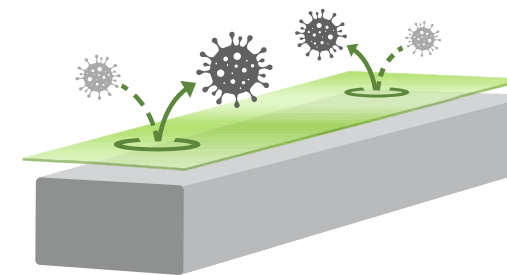

Bed components are individually powder-coated with a durable, anti-microbial polyester finish.

The bed adjusts vertically with virtually no horizontal motion, saving space and protecting the bed environment.





## SPECIFICATIONS

Height Floor to Mattress Support	7-7/8" – 30" / 200 mm – 762mm
Deck Angles – Back Elevation	0 – 70°
Deck Angles – Knee Elevation	0 – 30°
Bed Length – Overall Length	84-3/4", 88-3/4" / 2134 mm, 2235 mm
Bed Length – Sleep Surface	76", 80" / 1930 mm, 2032 mm
Bed Width – Sleep Surface	35" / 889 mm
Weight of Bed (without accessories)	160 lbs / 73 kg
Safe Working Load	500 lb / 227 kg

## WARRANTY

Electrical	5 years
Frame	15 years
Welds	15 years

## CERTIFICATIONS

Electrical & mechanical Safety Certification CSA C22.22 No. 601.1 (2005)
UL 60601-1 (2003)
CSA C22.22 No. 60601-2-38 (2003)
IEC 60601-2-38 (1996)
IEC 60601-1-2 (2004)

## ITEM CODE DESCRIPTION

QD1000	Advantage bed 35" x 80"
QD1000MLF	Advantage bed 35" x 80"; 360° mobility; FEC; floor lock
QD1000Y	Advantage bed 35" x 76"
QD1000YMLF	Advantage bed 35" x 76"; 360° mobility; FEC; floor lock
Q6599	Head end roller bumpers
RM930	Three position pivot assist rail - universal
RAD959	Molded Pivot Assist pair
RM950	Fixed assist rail - universal
RAD306H	Standard half head end side rail - set
RAD308H	Set of molded half head end side rails
RAD308H-L	Left side molded half head end side rail
RAD308H-R	Right side molded half head end side rail
RAD501	3/4 rail, 35" fixed width, folding, pair
Q6730	Trapeze Pole, Offset
Q6733	Mounting Frame, Trapeze
Q6587	Transport Dolly Set
CS82	IV pole adjustable
BS79B	I.V. Pole
RAD959	Pair Molded Pivot Assist Rails, quick width adjustable and frame mounted
Q6736	4" Extender Kit
Q6633	39" Fixed Expander Kit
Q6709	42" Fixed Expander Kit
Q6681	Battery Back-Up - factory Installed
Q6674	Foot end controller
Q6615	Floor Lock Kit



Get more information

800.268.8671(Span)

savariapatientcare.com



Span-America

| 800.268.8671

| PO Box 5231 Greenville SC 29606

Span-America is a division of Savaria Corporation, a global leader in accessibility. Learn more about Savaria at savaria.com

PL Span Advantage Bed CO# 2937 Rev. 0

©2023 Savaria Corporation. Savaria makes every effort to publish accurate information. Specifications and options are subject to change without notice.