


RHYTHM MULTI

Rhythm Multi Mattress System is a medium to high risk system for patients who are at risk of developing pressure ulcers according to your sound clinical judgment. Rhythm Multi can also be used for patients who have an existing stage 1, 2, 3 and 4 pressure ulcers, in conjunction with your hospital's policy on pressure area management.

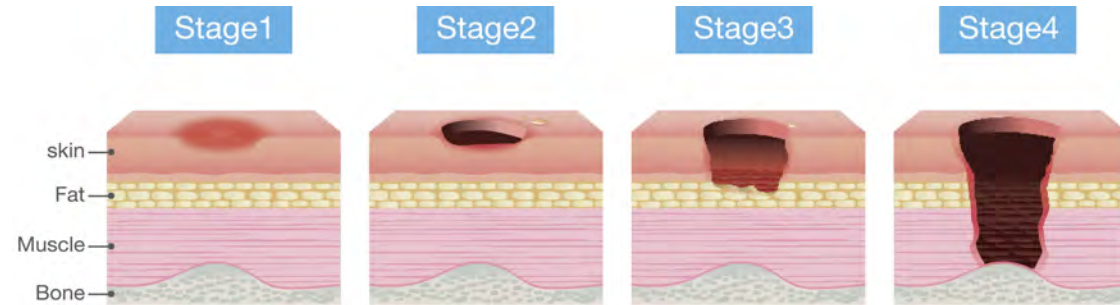

Optional Raised Bolster Cover available for added patient safety. Bolster cover has center access providing easy ingress and egress.

MASTER CONTROL UNIT	
Dimension:	(L)17.7"x (W)6.8"x (H)10.8"
Weight:	13.2 lbs
Rating Voltage:	110-120 volts
Max Current:	5A
Rating Frequency:	60Hz
Cycle Time:	3~95 mins 3 in 1 alternation
Min/Max Pressure:	10~61 mmHg +/- 3~4 mmHg
Protection Type:	Class I, Type BF Power outage visual and audio alarm Static / Pulsation / Alternating Pressure Auto-firm for fast mattress inflation Keypad lock-out function
Air flow rate:	> 1200 liters per minute

MATTRESS REPLACEMENT	
Model No:	Rhythm Multi
Dimensions:	(L)80"x (W)36"x (H)10"
Bariatric Size:	42, 48 inch width
Weight:	28.6 lbs
Weight Capacity:	600, 800, 1000 lbs
Cell Material:	Black Nylon with PU coating
Cell Number:	18 cells
Cover Material:	Blue Nylon with PU coating finish
Base Material:	Black Polyester with PVC coating CPR Quick Release System Cable Management System on bottom cover
HCPSC Code	E0277

PRIUS HEALTHCARE USA
Creating innovation and difference

T (813) 854-5464
E info@priushc.com
W priushcusa.com

160 SCARLET BLVD
OLDSMAR, FL 34677

RHYTHM MULTI

Provides you, the caregiver, with peace of mind, and your patient with world class care and comfort




PRIUS HEALTHCARE USA
Creating innovation and difference

priushcusa.com

Perfect for patients and caregivers




HAPPYHEEL™ provides extra comfort settings for the heels

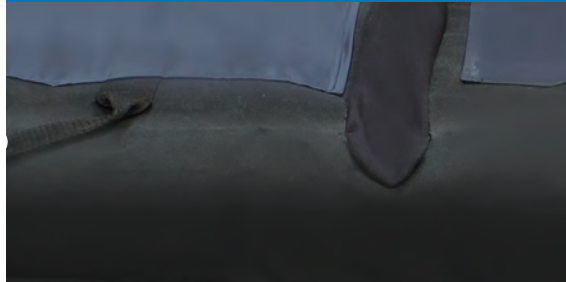

Low slip nylon under mattress base cover to reduce shearing and friction between mattress and bed frame



Nylon woven fabric with PU coating finish

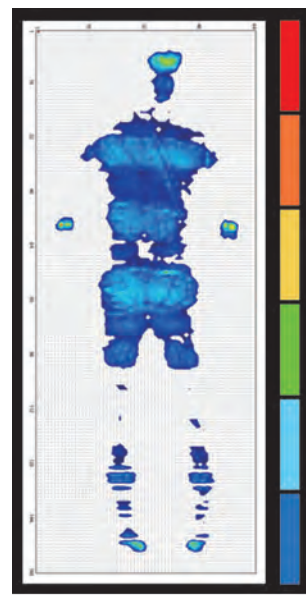


PRESSURE MAPPING

Rhythm Multi provides constant low pressure by redistributing the patient's pressure over a greater area than a standard mattress.



Foam Safety Mat to prevent bottoming out in case of power loss

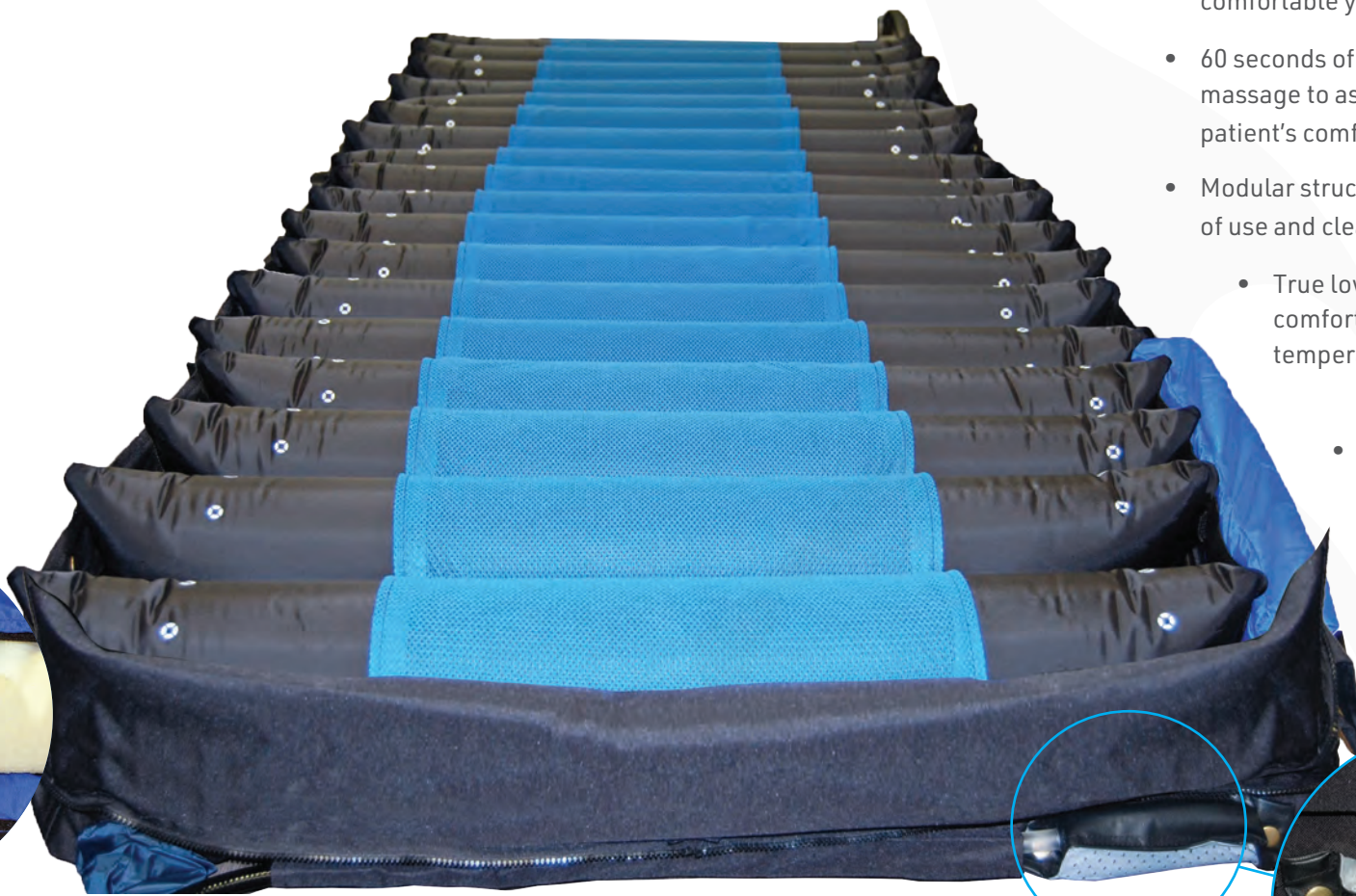

RHYTHM MULTI


- The Rhythm Multi replacement mattress can accommodate patient weights up to 1000 lbs
- 3-in-1 alternating pressure covers 66% of patient's body to provide a more comfortable yet still effective therapy
- 60 seconds of pulsation simulates a massage to assist in maximizing a patient's comfort
- Modular structure design ensures ease of use and cleaning
- True low air loss provides extra comfort for moisture and temperature control



- provides 3 selectable softness settings in the last 4 air cells to provide protection for patient's fragile heels



Exhaust port removes excess heat from pump lowering mattress surface temperature

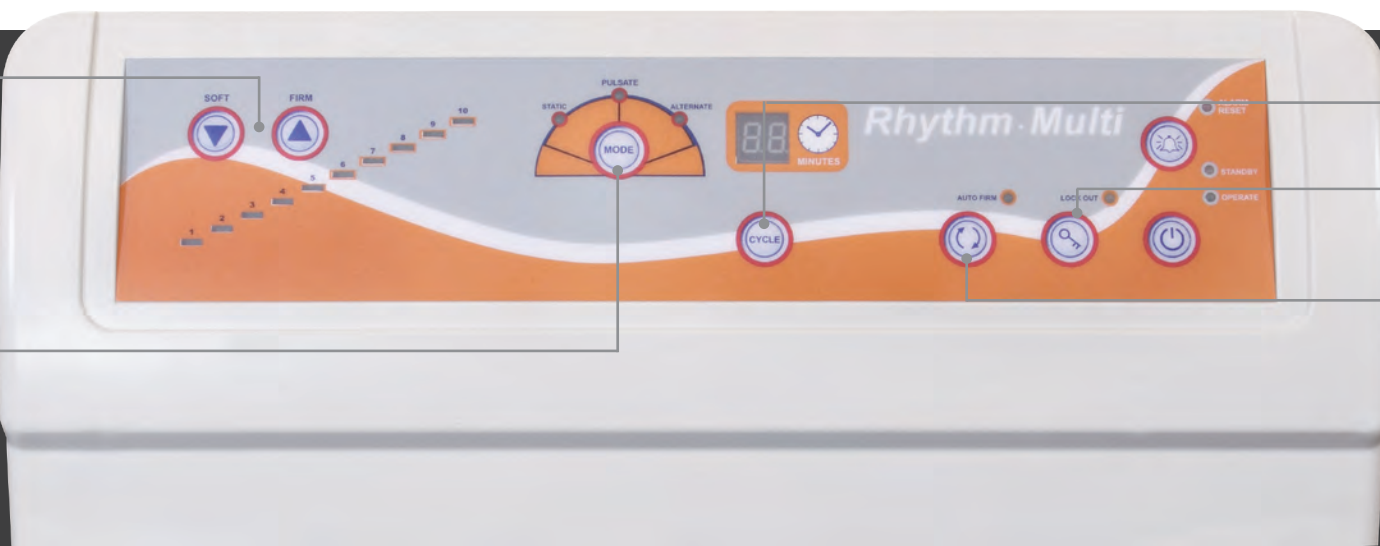

10 adjustable comfort settings provides different levels of comfort depending on the patient's therapy required

3 different modes of clinical therapy: static/pulsation/alternation

Alternation cycle times from 3-95 minutes

Control Panel Lock Out

Auto Firm/Patient Care control



Rhythm Multi is your partner in protection and treatment of pressure ulcers



Cable management keeps electrical cables off the floor



CPR quick release system deflates the system within 20 seconds



Hanging bracket fits on all type of beds