

Virtuoso Performance in Resident Care Beds










GlideAlign™ Retractable Deck Technology

Perfectly aligned with your facility's safety goals.

Superior fall prevention: By helping the user maintain vertical alignment, Encore helps ensure that all points in the bed environment remain more safely within user's reach (see illustrations p. 5).

Superior staff protection: Helps offset the tendency to slide down in bed as head of bed (HOB) is elevated, reducing the need for frequent repositioning that can lead to staff patient handling injuries.

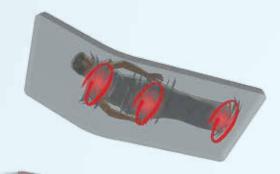


Encore's deck moves — so the user doesn't have to.

Superior protection from shearing

By reducing sliding during HOB elevation, helps minimize opportunity for skin shearing at vulnerable sacrum, scapula and heels.





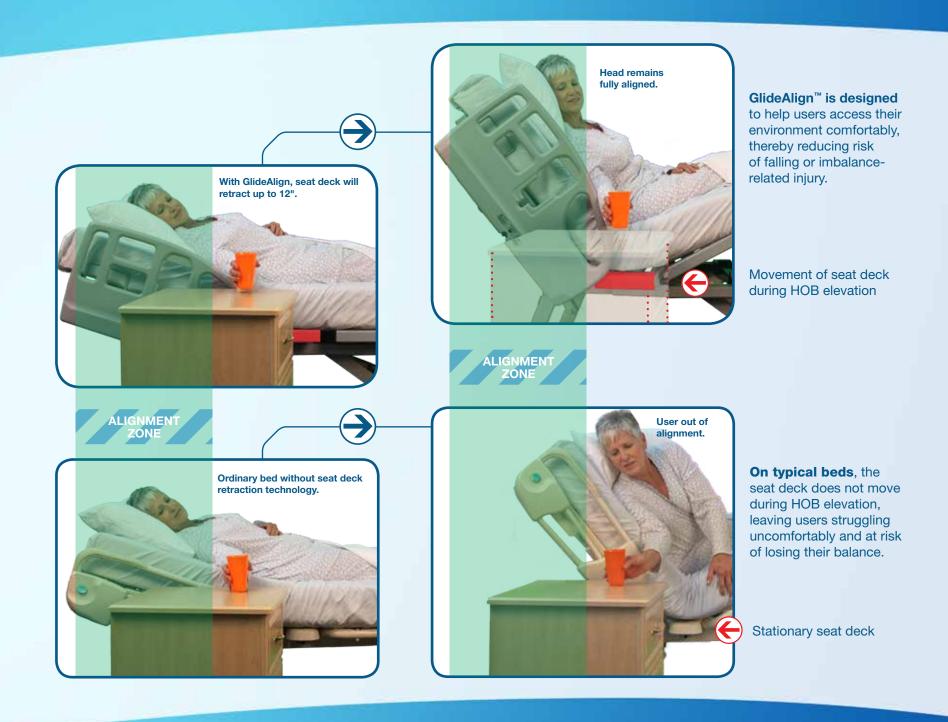

ReadyWide™ Integrated Expansion



No bolts. No tools. No pins. Retracts into deck. Nothing to store, lose, or become loose.



Tool-less, "two-click" expansion.



Navigates narrow doorways: Fast, hassle-free bed mobility.

- Simple, push-button operation: Click once to pull expanders out to 39" wide. Click again for 42" width.
- Can be factory-installed at time of order, or ordered separately as needed (see back page for options).
- Expands treatment flexibility, supported by Encore's 600 lb. Safe Working Load (SWL).



Wider surface facilitates enhanced user comfort, fall prevention, transfers and ADLs.





Smart Stop[™] low position auto-pause





Dual underbed indicators flash as bed enters "low height caution zone."

Improves low-bed safety and ergonomics:

- Automatically stops lowering when deck reaches 12.5" from floor.
- Lowering can then be resumed via foot-end (staff) control only.
- Provides cue for caregiver to confirm underbed clearance.
- Eliminates potential hazard of patient lowering bed onto unseen obstructions via hand pendant.
- Can aid in reducing injury, falls by users who might otherwise attempt egress from a "too-low" position.

90° angle

Helps ensure safer, more ergonomically correct transfers. Positions typical user

at ideal height* and leg angle for egress/regress. (*12.5", plus height of mattress)





Key Standard Features

Head/footboard options see page 10

Mattress Retainers

Durable flexible polyester finish on individually powder-coated parts.

Four (4)
synchronized DC
self-lubricating
actuators





Precision welds with lifetime guarantee.

Baseboard bumper

Hospital-grade ball bearing swivel casters, provide full mobility at any height. Head-end includes directional lock-out.



Dual-position, dualintensity underbed lighting.

- Standard on both left and right.
- Selectable bright or discreet nighttime floor level illumination.



Travels virtually straight up and down; maintains ample underbed clearance.







Backlit 9-button pendant with side-to-side adaptability.

CPR Button

Advanced Clinical Positioning



Reverse Trendelenburg



Vascular foot ratchet elevates feet.

Standard model can be used at 76" or 80"L.



Auto Contour



Comfort Chair



Trendelenburg

 Smart Stop[™] low height warning/override



Foot-end floor lock

can be custom-ordered for head end.



Options & Accessories

Finishes



Standard Profiles

Essex Lincoln Kent Standard



Elite Profiles

Lethbridge Pembrooke Trenton Princeton



Recessed pendant holder,

included, fits all rails and supports pendant facing either toward user or toward caregiver.



Bed extenders to 84"

Q6690






Transport Dollies for easy mobility and passage through narrow doorways; simplified access to under-bed components.

Q6699 (1 set)







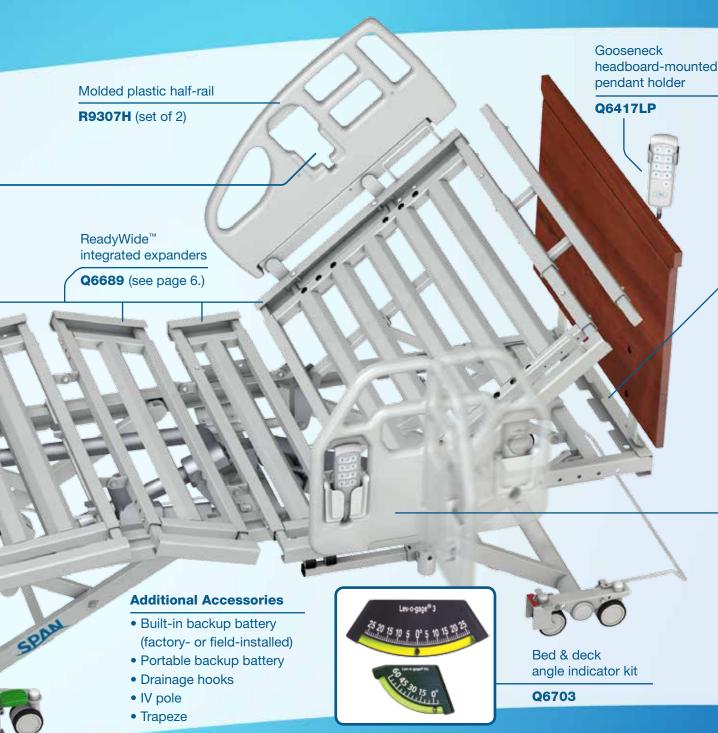

Honeycomb roller bumper kit

Q6708



Ergonomic 3-position assist device locks perpendicularly for sturdy ingress/ egress support, protecting both user and caregiver. Both 3-position and half-length assist devices are made of molded, pleasing-to-touch plastic.

R9955 (1 ea, interchangeable L or R) **R9956** (set of 2)

	Alighted			in the second of				
MODEL NUMBER	Glider and	Small dual u	600th ookon	Floot or lo	600 Knorth	Read	Built's institution	
MC9R	X	X	X X	X	X			
MC9RW	Х	X	X	Х	Х	Х		
MC9RZ	Х	Х	Х	Х	Х		Х	
MC9RWZ	Х	Х	Х	х	х	Х	Х	

	(a) - a
Height Floor to Mattress Support	7-7/8" - 31"
Deck Angles - Back Elevation	0 – 70°
Deck Angles - Knee Elevation	0 – 30°
Bed Length - Overall Length	82", 86", 90"
Bed Length - Retraction Travel	12"
Bed Length - Sleep Surface	76", 80", 84"
Bed Width - With Side Rails	40-1/2", 43-1/2", 46-1/2"
Bed Width - Sleep Surface	35", 39", 42"
Mattress Support Deck Const.	Welded steel tube and slats
Mattress Stops	Steel Corner and Side Stops
Bed End Panels	Permanent or removable mounting
Electrical & Mechanical Certification	ANSI/AAMI ES60601-1:2005, Medical electrical equipment
	CAN/CSA-C22.2 NO. 60601-1-08 - Medical electrical equipment
	CAN/CSA-C22.2 NO. 60601-2-52:11 - Medical Beds
	CAN/CSA-C22.2 NO. 60601-1-2-08 - EMC
Classification	Class 1, Type B
Classification VAC/Hertz/Phase	Class 1, Type B 120/60/1

Cord Type/Length	2.5 Meters minimum, 3 prong hospital grade plug		
Electric Motors	Four (4) DC self-lubricating linear actuators		
Operation	10% Maximum Duty Cycle (2 minutes on,18 minutes off)		
Emergency Power back-up battery Pack	24V DC Internal or External Battery with Charger		
Ingress Protection Rating	IPX4 (Electrical components IP66)		
Audible Acoustic Energy	Less than 80 dBA		
Weight of Bed	253 lbs / 115 kg without any accessories		
Safe Working Load	600 lbs / 272 kg (Includes weight of accessories, options and bedding)		
Mobility	At any height (7-7/8" - 31")		
Swivel Casters	Heavy duty dual 3" Casters		
Warranty	Mechanical: 5 years; Frame: 10 years; Welds: 15 years.		

Outfit your Encore with one of Span's many outstanding support surfaces up to 42" wide. Options include:

PressureGuard Easy Air®



Geo-Mattress Atlas® Bariatric Foam



Geo-Mattress® UltraMax™



Custom Care® Convertible



